## Bionic Navigator™ 3D System Software Upgrade

Installer for BN3.1 (SC-7102)  
Installation Guide

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90970819-02 REV C

### **Guarantees**

Boston Scientific Corporation reserves the right to modify, without prior notice, information relating to its products in order to improve their reliability or operating capacity.

Drawings are for illustration purposes only.

### **Trademarks**

All trademarks are the property of their respective holders.

### **For Clinician Programmer Models:**

<b>Model Number</b>	<b>Description</b>
NM-7153-30	Clinician Programmer
NM-7153-30R	Clinician Programmer (Refurbished)

# Installer for BN3.1 (SC-7102) Installation Guide

## Description

The enclosed Flash Drive contains an executable file that will upgrade the user's Bionic Navigator™ 3D Clinician Programmer software version to 9028372-110.

## Package Contents

(1) Installer for BN3.1 containing installation file for 9028372-110

## Software Installation

### Boston Scientific Neuromodulation (BSN) employees:

To complete the installation, you will need your Clinician Programmer (CP) and another VPN-enabled BSN computer. This upgrade may only be performed on a compatible model CP, listed in this manual. On your CP, installing this update will upgrade your Bionic Navigator software to 9028372-110. On your BSN computer, after completing the VPN process, you will follow steps outlined below to 1) receive an activation key and 2) provide a confirmation key to allow BSN to track the upgrade. The steps below must be repeated to install 9028372-110 on each compatible CP as each installation requires a unique activation key.

If you encounter any issues during installation, please call Patient Care at 866-360-4747.

### Non-BSN employees:

To complete the installation, you will need your Clinician Programmer (CP). This upgrade may only be performed on a compatible model CP, listed in this manual. On your CP, installing this update will upgrade your Bionic Navigator software to 9028372-110. To obtain activation keys, provide confirmation keys to BSN, and to report any issues, please contact your local sales representative. The steps below must be repeated to install 9028372-110 on each compatible CP as each installation requires a unique activation key.

en


### Software Installation on CP

Verify that you have received the 9028372-110 Flash Drive by following the steps below:

1. Open pouch containing the Flash Drive.
2. Turn on CP.
3. Navigate to CP desktop and ensure that all applications are closed.
4. Insert Flash Drive into an unused USB port.


5. Navigate to the installation file on the Flash Drive.

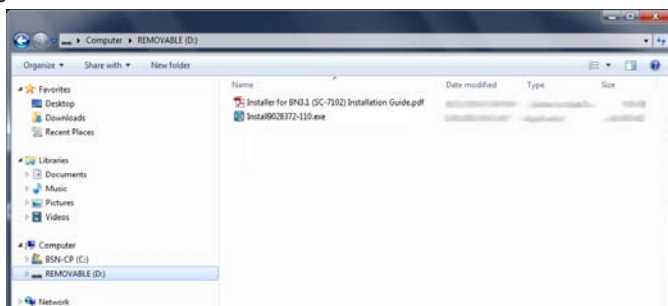6. Verify that the file displayed reads "Insta9028372-110" as shown above.

7. Open the installation file and click Next on the launch screen to start the installation process.

en

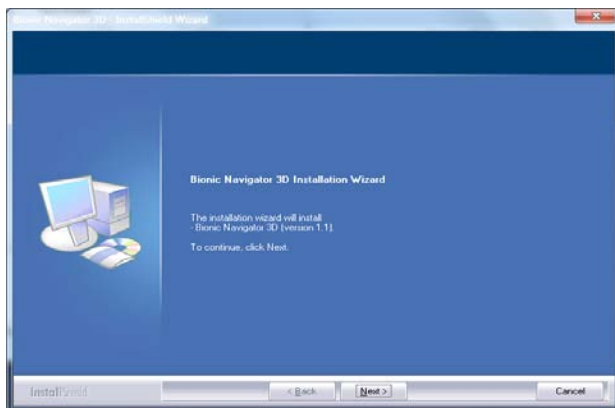8. Verify that the displayed Software Part Number is 9028372-110, as shown below:

**Note:** Obtaining an activation key is described in the following steps.


*BSN employees: continue to steps 9-14.*

*Non-BSN employees: Contact your local sales representative to obtain an activation key, then continue to step 14.*


en

## Web Activation Portal

9. On a BSN computer, log onto the BSN VPN.
10. Once connected to the VPN, open a web browser and connect to the BSN Product Activation Portal using the following URL:  
<http://valapp16/SAP/ActivateSAP/>

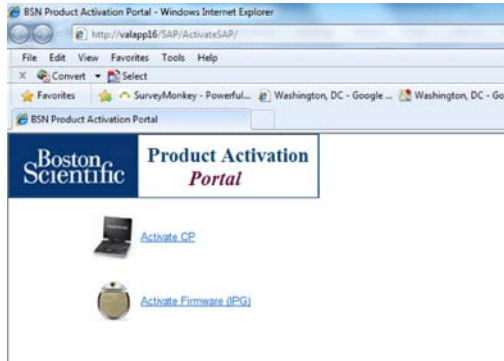11. Click on "Activate CP".

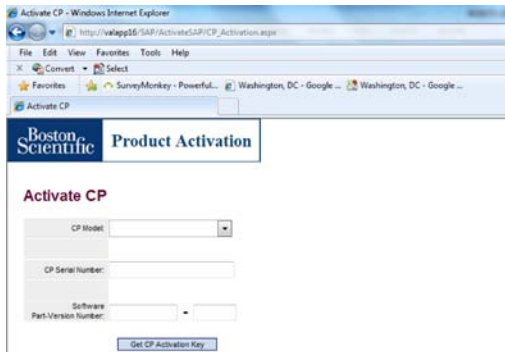
12. Use the drop down arrow to select the CP Model number, then input the CP Serial Number and Software Part Version Number displayed on the CP screen into the appropriate field on the web page.

**Note:** Serial numbers are case-sensitive. Note also that the letter "O" and the number "0" look similar and must be distinguished.

13. Click “Get CP Activation Key” to generate the activation key. Record the activation key, to be used in the following “Software Installation on CP” step. The activation key is valid for five (5) days.
14. Click on the “Return to Product Activation Portal” link.


## Software Installation on CP

15. On the CP, input the activation key generated from the web activation portal.
16. Click the Activate button. This process may take several minutes as the system is being reconfigured to support the the new application. Various status messages will be displayed as the Installer updates, including the following:

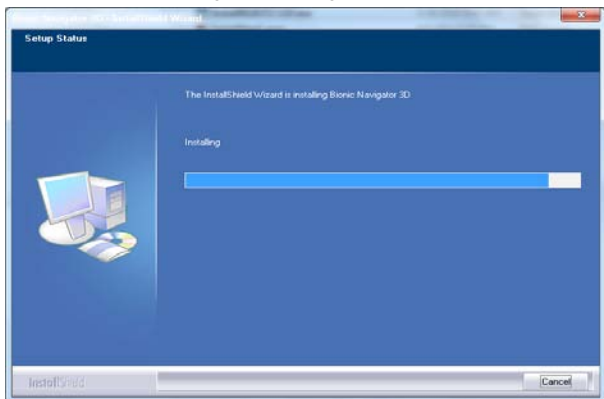en

17. After successful installation, a Confirmation Key will be displayed.

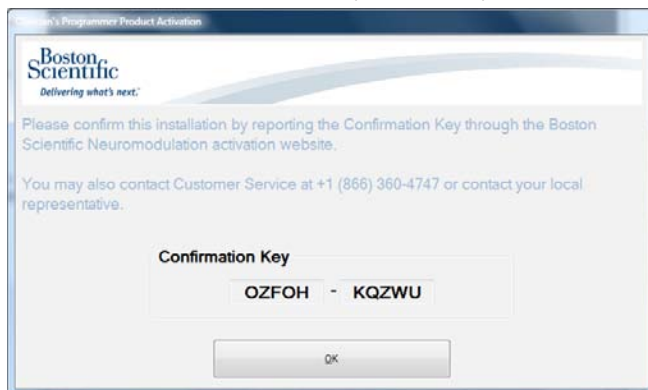18. Record the confirmation key displayed on the screen.

**Note:** The Confirmation Key is also displayed on the Configuration->Software Information tab of the Bionic Navigator 3D system, as shown below.

BSN employees: continue to steps 19-23.

Non-BSN employees: Contact your local sales representative to provide the Confirmation Key. Installation is now complete, and the Flash Drive may be removed from the USB port.

19. Open Bionic Navigator 3D and navigate to Configurations Tab.


20. Click on the Software Information Tab [computer icon] and view the Software Part Number. It should now read: "9028372-110", as shown in the following screen:

**Note:** If you do not see a successful installation message, or the Software Part Number does not display correctly, contact your local field support representative.

## Web Activation on BSN Portal


21. On the Product Activation Portal, click on "Confirm Activation".

22. Input the confirmation key from the CP into the displayed field on the web activation portal.

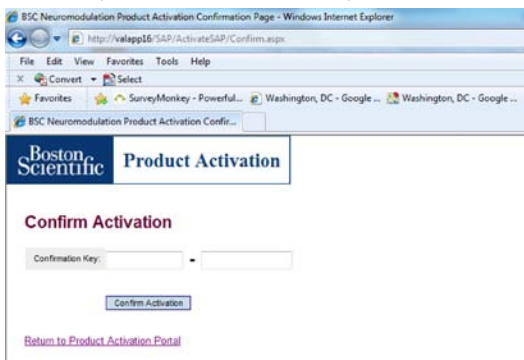en

23. Click “Confirm Activation” button.
24. Click Finish to complete installation.

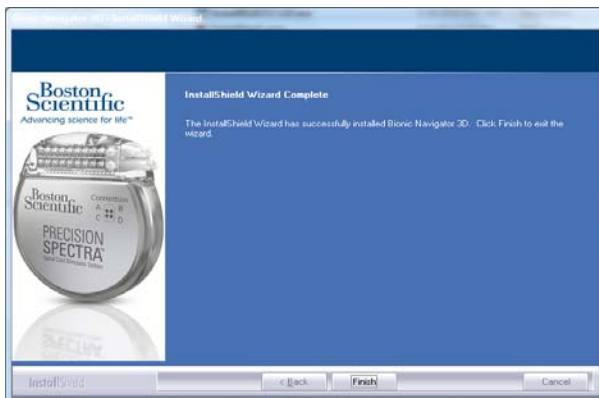The installation is now complete. The Flash Drive may be removed from the CP.

---

## Customer Service - USA

In USA, to order replacement parts, or for answers to any questions you might have, please contact the Customer Service Department:

- Phone: (866) 360-4747,  
+1 (661) 949-4747
- Fax: (661) 362-1503

Boston Scientific Neuromodulation  
25155 Rye Canyon Loop  
Valencia, CA 91355, USA

Outside USA, if you have a specific question or issue, please contact your healthcare professional.

## Customer Service - Outside USA

If you have any other questions or need to contact Boston Scientific, choose your locality from the contact list at the end of this document.

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or its affiliates. All rights reserved.
Boston Scientific (Australia) Pty Ltd
Boston Scientific International S.A.
© 2015 Boston Scientific Corporation